

## ENDOWMENT

- 1x — Unit LogiSTIM® EM-19W
- 2x — Wires for electrodes connection
- 4x — Self-adhering electrodes 50x50 mm with connector 2 mm socket type.
- 2x — Batteries type AA1.5 V (LR6)
- 1x — Instructions booklet
- 1x — Bag with compartments

Optional:  
kit for visage treatments



[www.beac.it](http://www.beac.it)

## WARNINGS!

It is a medical device.  
Carefully read the warnings and the instructions for use. The unit cannot be used: - on patients with cardiac pacemakers; - on pregnant women - on areas affected by inflammatory processes.

## TECHNICAL SPECIFICATIONS

**Supply:** 2 alkaline batteries 1,5 V type AA (LR6);

**Output:** 100 mApp over 1000 Ohm load (using 200  $\mu$ S pulse-width)

**Pulse width:** From 50 $\mu$ S to 400 $\mu$ S adjustable in steps of 50  $\mu$ S

**Frequency:** from 1Hz to 110Hz (1-2-3-4-5-10-12-14-16-18-20-25-30-35-40-45-50-60-70-80-90-100-110).

**Wave shape:** Bi-phasic symmetrical pulses – Bi-phasic alternated pulses.

**Treatment timer:** Continuous, 10min, 20 min, 30min, 45min, 60min, 90min.


**ACTION/REST time:** From 1 to 40 sec, adjustable in steps of 1 sec.

**RAMP time:** From 1 to 5 sec, adjustable in steps of 1 sec.

**Dimensions:** 100mm x 68mm x 35mm.

**Weight:** 150 g. including batteries

# LogiSTIM® EM-19W

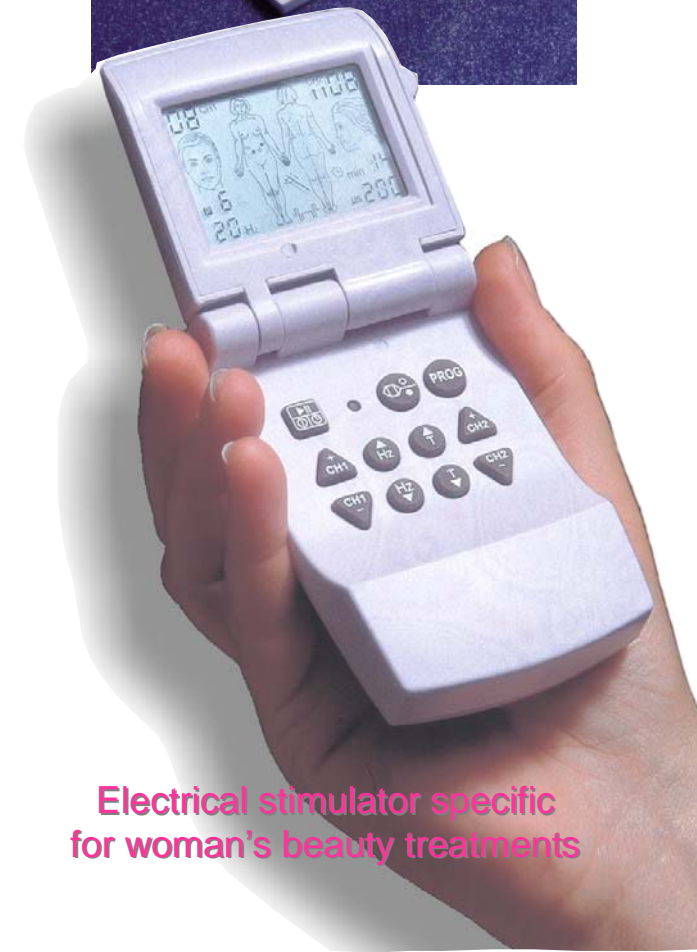

Electrical stimulator specific  
for woman's beauty treatments

Distributed by: **BEAC**  
27040 – Portalbera (PAVIA) ITALY  
Certified Quality System - EN ISO 9001 – ISO 13485



*“tecnoLogic beauty products, made for the modern woman”*

## Just select the area, then entrust to...

# LogiSTIM® W

### How does LogiSTIM® W work?

Through a couple of self-adhering electrodes, stimulation pulses are applied over the muscle; as result the muscle will alternately contract and relax, exactly as it happens during the physiological activity but with the advantage of avoiding the general sense of fatigue. Given the importance of the correct electrodes positioning, LogiSTIM® offers a UNIQUE feature on the market; to ensure the best effectiveness of the treatment and simplify the use, the right position of any electrodes is automatically shown on the display, according to the selected area to be treated. The user is only asked to select the area, then increase the stimulation level and progressively increase the frequency for improving the comfort. By selecting one of the 16 programs in memory, the following parameters will be automatically set:

- Treatment area,
- Right electrodes position,
- Wave shape,
- Pulse width.

**Personal electrostimulation unit, specifically designed for beauty/fitness treatments (body/visage), with easy graphical selection of the area.**



**Direct selection of the area to be treated**

**Best electrodes positioning, automatically shown on the display**

**Function IntelliSTIM® with automatic interactive parameters adjustment**

**Any parameters shown, adjustable in real time, during the session.**

**Automatic selection of the wave-shape, according to the area**

**4 free programs**

- Lift-UP the BREAST (m. pectoral)
- Slim ABDOMEN (m. upper abdominal)
- Trim the WAIST (m. abdominal + oblique)
- Lift-UP lower GLUTEI
- Re-define upper GLUTEI
- Improve medium GLUTEI
- Reduce THIGHS
- Re-style ARMS
- Reinforce the BACK
- Tone-up HAMSTRINGS
- Improve LEG circulation
- Enhance the CALVES
- Reduce MOUTH WRINKLES
- Reduce EYES WRINKLES
- Reduce DOUBLE CHIN
- Lift SAGGING CHEEKS