#### **PACKAGE**

1x — Unit BE-28T

2x — Wires for electrodes connection

4x — Self-adhering electrodes, with 2 mm socket connection

2x — Alkaline batteries AA1.5 V (LR6)

1x — Instructions booklet

1x — Stuffed bag with compartments



### WARNINGS!

BE-28T is a medical device CE.

Rev. 1 of 13/08/2006

Carefully read the warnings and the instructions for use. The unit cannot be used: - on patients with cardiac pacemakers; - on pregnant women - on areas affected by inflammatory processes.

#### **TECHNICAL SPECIFICATIONS**

Supply: 2 alkaline batteries 1,5 V AA (LR6)

Maximum 130 Output: mApp across 500 Ohm load

Pulse Width: From 50µS to 250µS adjustable in steps of 10µS

Pulse rate: From 1Hz to 150Hz (1-2-3-4-5-10-12-14-16-18-20-25-30-35-40-45-50-60-70-80-90-100-110-120-130-140-150)


Waveform: Symmetrical Bi-phasic rectangular - Asymmetrical Biphasic rectangular - Mono-phasic rectangular

Treatment Timer: Continuous. 15min, 30min, 45min, 60min, 90min.

Dimensions: 138mm x 6 8mm x

28mm

Weigth: 160 g. including batteries

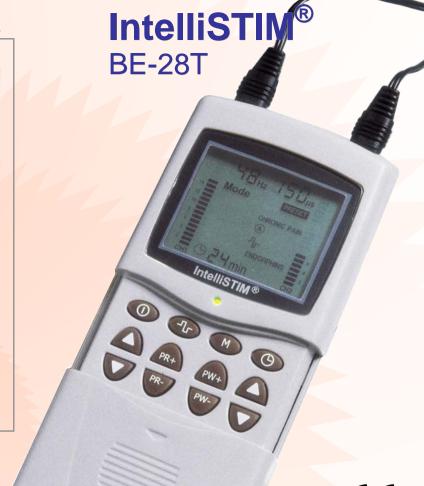

Professional unit for pain management through electrical stimulation (T.E.N.S.).

Distributed by:

**BEAC** 

27040 - Portalbera (PAVIA) ITALY **CERTIFIED Quality Management System** 



ISO 9001 - ISO 13485

# We exactly know what you expect by a professional analgesic unit...

## ...and we developed...

By a modern analgesic electrical stimulation unit (T.E.N.S.) a professional user expects that it is: **effective** (for achieving therapeutic aims), **easy to operate** (for avoiding the user to waste his time), **flexible** (for allowing to treat any pain situation), **robust** (to resist for long time). Further to the above, **IntelliSTIM® BE-28T** also offers:

- > Large LCD display
- > Intuitive Operation
- > 9 operating programs (editable)
- > any parameters shown
- > "edit & feel" function
- > ACUTE and CHRONIC pain protocols
- > 3 Selectable wave shapes





# Professional unit for pain management.

- Dual intensity-independent channels
- Last used program loading
- Count-down timer with switch-off at the end
- Auto-switch-off if not in use
- Accumulative recording of the treatment time (in hours)
- Electrodes protection circuit (just in the model BE-28TC)
- Lock function

# ENDORPHINS or "Gate -Control"?

Among the other parameters, on the display, it is also shown the applied analgesic technique, according to the set-up of the unit.